



Illuminator aka Inner-Luminator



Low-Level Laser Therapy accelerates the body's natural healing process through photo-bio-modulation.



Clinical Benefits:

- Non-pharmacologic pain relief
- Effective for difficult conditions
- Alternative to surgery
- Fast treatment times
- Easy, comfortable, non-invasive treatment
- Scientifically-proven

Practice Benefits:

- Provides a treatment modality that will help you attract more referrals
- Provides an alternative or supplement to prescription drugs & surgery
- Increase revenue with this versatile modality

Inner-Luminator—uses low-level light energy to treat a wide array of conditions, ranging from:

temporomandibular disorder (chronic jaw pain)
sprains
tendinitis
back pain
wound healing

arthritis
carpal tunnel syndrome,
diabetic neuropathy
fibromyalgia,
leg ulcers



Statistical probabilities:

The UK, USA, North America, Australia, and New Zealand, have the highest recorded rates of breast cancer in the world. While most cases are found in women over 40, 25% of all cases are women under 50, and the incidence of breast cancer in much younger women seems to be increasing.

Used by some physical therapists, physiatrists, chiropractors, and occupational therapists, it is different from conventional "hot" laser surgery.

Inner-Luminator has anti-inflammatory effects, which may help repair connective and other tissues and release pain-relieving endorphins.

Breast Health


Inner-Luminator is a health and wellbeing product for women that can give a rapid, accurate, scientific, and non-destructive indication of breast disease, especially for carcinoma, breast cancer.

How does it work?

Using the latest medical technology of Trans-Dermal Illumination,

Inner-Luminator works by shining a powerful light through your breast.

This light is absorbed by Hemoglobin and makes the veins in your breast appear as dark lines. Tumors and malignant lesions will show up as dark areas due to the blood vessels that form to feed such malignancies. Harmless fibroid cysts will allow the light to flow through them, but potential abnormalities will block the light beam, revealing dark spots. Such dark areas around malignancies should show before you can feel any lump or abnormality by hand self-examination, thus giving you an early warning to consult your doctor.

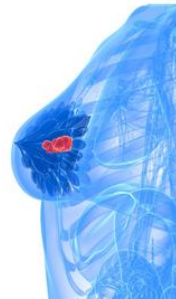

Peace of mind

There are no harmful side effects from using **Inner-Luminator**, so you can use it regularly for breast self-examination on a precautionary basis. Should any abnormality show you would have the benefit of very early warning, which greatly enhances the chances of successful treatment. You do not have to remember what previous examinations looked like, just look for changes that show as dark areas.

Many lumps found in breasts can be harmless, and benign breast tumors do not spread outside the breast but may increase the chance of developing cancer.

Inner-Luminator in particular will show abnormal blood vessel changes at an early stage.

Other Breast Cancer awareness points to look for in self-examination.

If you notice any of these changes, or dark areas by using a **Inner-Luminator**, you should consult your doctor for further evaluation.

Lumps

Change in the shape or size of your breast

Change in the skin, particularly dimpling or puckering

Breast pain that does not go away after a period

Change in the appearance of the nipple, or discharge from the nipple

Lumpy areas or thickening of the breast tissue

These statements have not been approved by the FDA. Results may vary. The device does not diagnose, prevent, treat or cure disease. Consult health care provider as necessary.

For purchasing, <https://shop.solexnation.com/JocelyneAPRN> or call (860) 830-1115



Illuminator

By Solex



illuminator™



The Solex Illuminator is an advanced topical wellness device designed with three specialized laser diodes.

BENEFITS:

- Stimulates collagen production for firmer, smoother, more radiant skin
- Deep infrared laser tissue penetration for
 - Pain Support
 - Vasal Dilation
 - Blood Cell Function
- Rejuvenates skin, and visibly reduces blemishes and other signs of aging*

*These statements have not been evaluated by the FDA. This product is not intended to diagnose, treat, cure or prevent any disease.

SOLEX 

SOLEXLLC.COM | 385.203.0033
SUPPORT@SOLEXLLC.COM

- **Infrared red-light therapy:** Enhances tissue rejuvenation and stimulates collagen to improve the appearance and vibrancy of your skin.*
- **Therapeutic warmth:** Increases the absorption of skin care products by enhancing skin permeability and blood flow.*

HOW TO USE:

- Put device against skin & move slowly in an upward direction
- Use in two minute increments in any one area
- Protect eyes from exposure to laser diodes

INTENSITY LASER SETTINGS:

- Intensity 1: 20mw • Intensity 3: 120mw
- Intensity 2: 50mw • Intensity 4: 180mw

CLEANING INSTRUCTIONS:


- Use an alcohol wipe or any antiseptic wipe to disinfect (Do not submerge in water).
- Air-dry the device or use a clean towel to dry.
- Make sure the device is completely dry before using it again.

SOLEX 

SOLEXLLC.COM | 385.203.0033
SUPPORT@SOLEXLLC.COM


These statements have not been evaluated by the FDA.

The Illuminator is not intended to diagnose, prevent, treat or cure any disease. Refer to Solex product use guidelines.

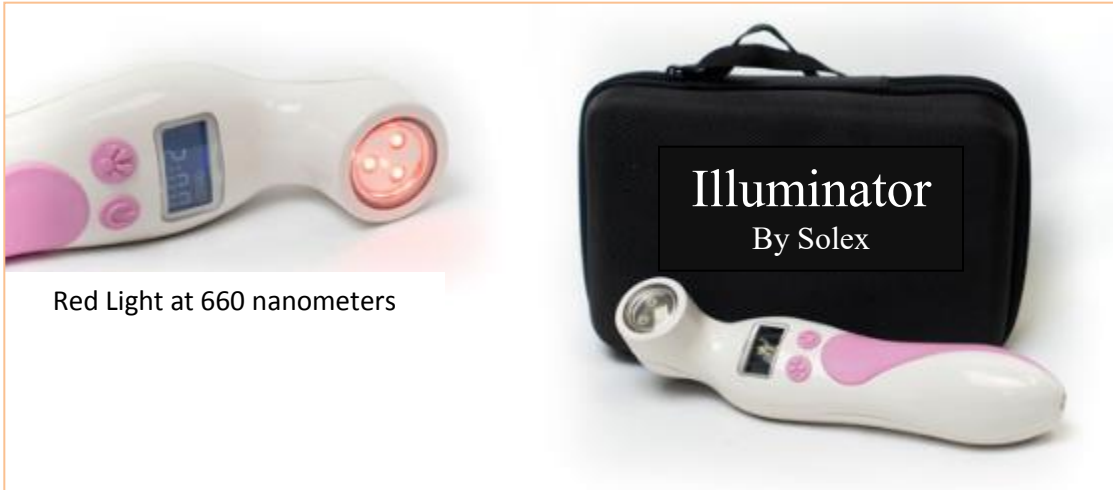

Red Light at 660 nanometers

The table below referenced from the book “*Red Light Therapy*” shows Conditions affected by specific wavelengths.
The Illuminator emits red light at 660 nanometers.

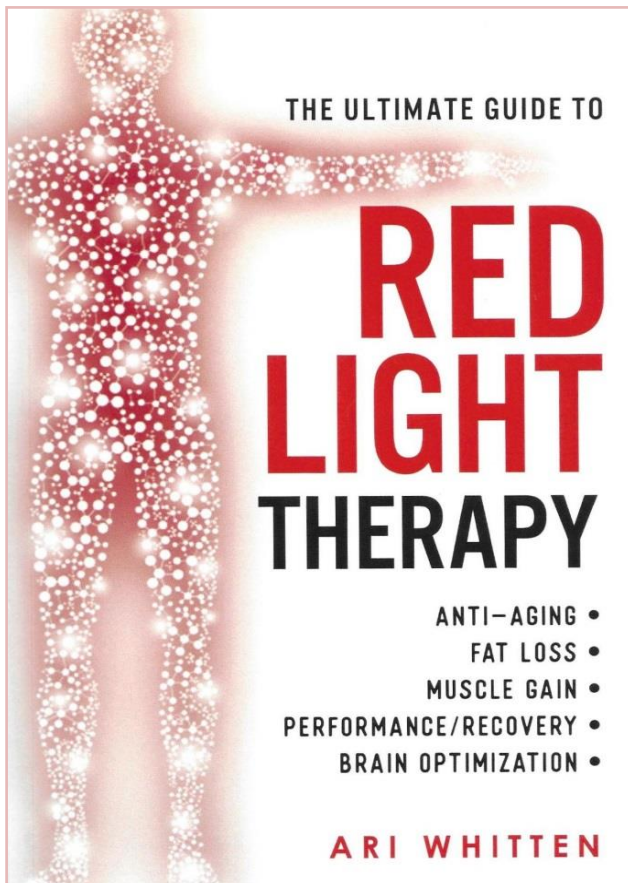

Condition	Wavelengths shown to benefit condition (in nanometers)
Alzheimer's	627, 670, 800, 810, 1070
Cognitive performance	660, 810, 1064
Depression	630, 650, 810
High blood pressure	660, 780
Eye health	630, 650, 660, 670, 680, 780, 810
Hair growth	630, 650, 830
Exercise performance	630, 640, 660, 670, 810, 850
Arthritis	630, 660, 785, 810, 830, 910
Oral health	630, 640, 650, 660, 670, 685, 780, 790, 810, 830, 850, 900, 940, 980, 1060
Pain	650, 660, 780, 810, 820, 830, 900, 980
Fertility	630, 650, 830
Skin rejuvenation	630, 650, 660, 830, 890
Thyroid health	630, 780, 830

These statements have not been approved by the FDA. Results may vary. The device does not diagnose, prevent, treat or cure disease. Consult health care provider as necessary.

For purchasing, <https://shop.solexnation.com/JocelyneAPRN> or call (860) 830-1115